



INTENSIVE CARE VENTILATION

elisa 300

Homecare
 Pneumology
 Neonatology
 Anaesthesia
 INTENSIVE CARE VENTILATION
 Sleep Diagnostics
 Service
 Patient Support

Intensive care ventilation.
Simple, effective
and lung protective.

Contents

elisa 300 The new compact range in intensive care ventilation with the latest turbine technology	4
The elisa principle – agile system design for individualised ventilation therapy	5
Instant View Technology in control at all times	6
Easy Access Bar precise operation even in stressful situations	10
A clean affair a simple way to prevent nosocomial pneumonia	9
PEEPfinder® gold standard bedside lung diagnostics	11
Tools to assist the weaning process there are no simple answers when weaning fails	12
Options & choices our modular system at a glance	14

elisa 300

The new compact range in intensive care ventilation
with the latest turbine technology



elisa 300



With a peak flow of up to 300 litres per minute, the high-performance, noise-optimised turbine guarantees sufficient flow capacities.

elisa 300 combines compact-class advantages with the performance characteristics of a modern universal ventilator for invasive and non-invasive ventilation therapy.

The innovative user interface and the comprehensive device configuration options are the basis for versatile application options in intensive care, intermediate care, emergency rooms or during in-hospital transport. The 12.1-inch display with a stunning colour performance is the key operating element to guarantee simple operation. Numerous functions provide support with daily routine tasks.

The elisa principle – agile system design for individualised ventilation therapy




Innovative. Intuitive. Sustainable.
The elisa family.

In theory, things couldn't be simpler - air has to go in and out.

Clinical practice, in contrast, consists of a wide range of requirements for modes, setting parameters, diagnostic options, and therapeutic manoeuvres. As a result, ventilators quickly become overly complex, subject to compromise or require intensive training. The agile system architecture of the elisa family implements the respective hospital standard to reduce operator errors and work time.

Whether customised rounds view, resuscitation mode or Fastwean support – the user interface can be perfectly configured to suit your needs.

The platform concept enables situation-based configuration. The flexible system architecture allows for integrating future requirements as well as medical and technical developments.

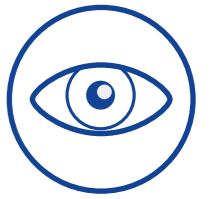

Instant View Technology

in control at all times

Don't miss the wood for the trees – instantly assess the current ventilation situation and identify developing problems.


Time is a scarce resource in everyday clinical practice. Increasing workloads, critical situations as well as normal routine place high demands on medical personnel. The cumbersome operation of complicated devices causes additional stress and creates sources of error. That calls for innovative technology which offers a clear overview of the required information in a structured format. Simply smart!

Instant View Technology

The Instant View Technology gives you an intuitive grasp of the patient's situation. Trends and necessary interventions are immediately evident. Deviations are clearly obvious without the need to read individual measuring values.



0% MV spontaneous 100%		
OFF	PEAK 16 mbar	40
	Pplateau 16 mbar	OFF
3	PEEP 5 mbar	7
	Leakage 3 %	50
6	RR 16 /min	18
400	VTe 325 ml	OFF
2.5	MVe 6.9 L	12.0
18	O2 30 %	26
	etCO2 mmHg	60
OFF	-	

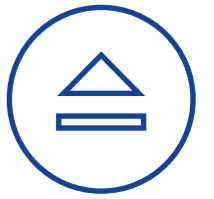

A clean affair

a simple way to prevent nosocomial pneumonia


The device's numerous individual functions and architecture support compliance with suitable infection prevention measures.

Pneumonia is the most common nosocomial infection occurring in ventilation patients. It leads to extended hospital stays and increases lethality by up to 30%.

The elisa 300 features a number of functions to support the necessary measures for reducing nosocomial infections. The design of the modern intensive care ventilators eliminates hygienic problem zones such as dirt-collecting corners or rotary knobs and allows for easy cleaning and disinfection. The Valve Bar comprises all elements that can be directly or indirectly contaminated via the respiratory tract and makes it easy to quickly replace all patient-side connections to effectively prevent cross-contamination.

The materials used guarantee continued functionality, even under the most severe conditions such as mechanical strain and repeated autoclaving.


The configurable hygiene function supports the implementation of internal hospital hygiene standards without the need for complex RFID technology or the purchase of expensive special tube systems. It comprises all potentially critical parts such as nebulizers, HME filters, tube extensions, and suction systems.

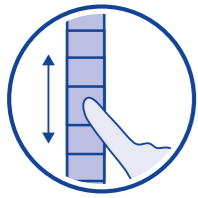

Easy Access Bar


precise operation even in stressful situations

New answers are needed for intelligent operation – the Easy Access Bar enables fast response.

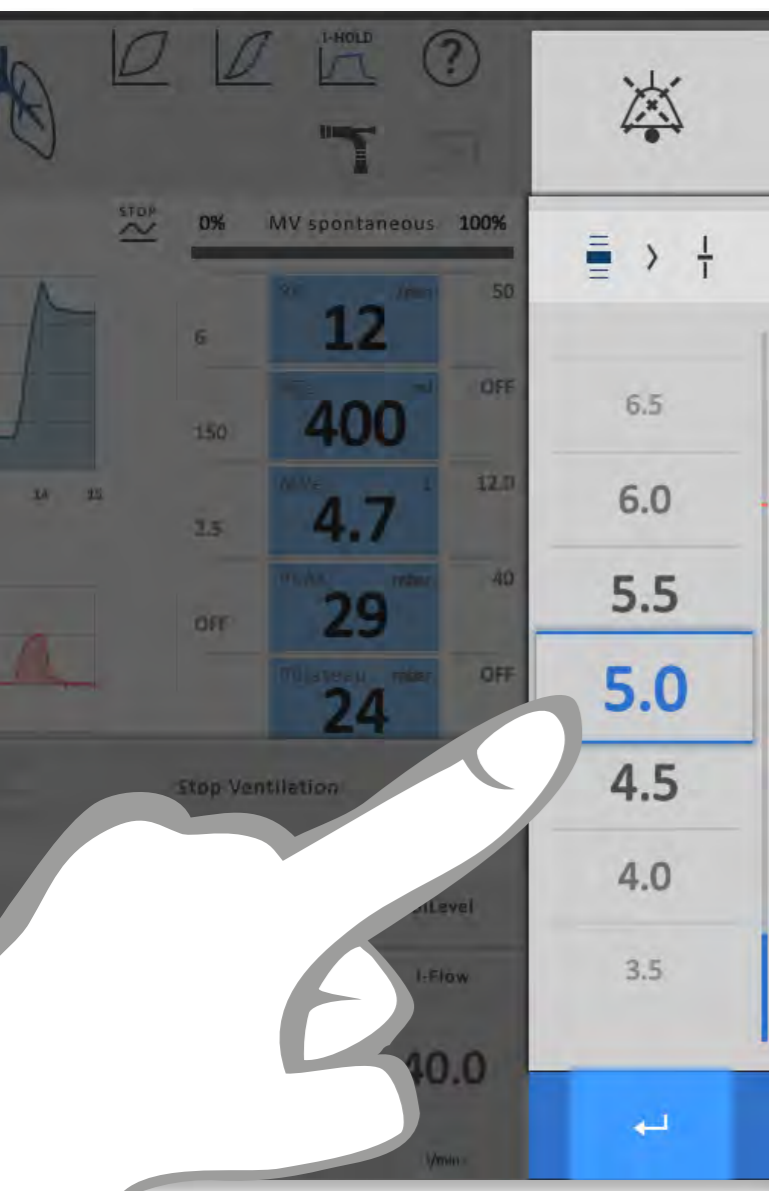

The Easy Access Bar of the intensive care ventilator elisa 300 lets you choose the required settings with precision and ease, even in stressful situations. The touchscreen operation provides intuitively understandable, unmistakable feedback on the selected setting. Since all numerical values and setting parameters are consistently arranged in the same location, operating the devices becomes an easy routine that does not fail in critical situations.

The absence of conventional rotary knobs makes operation easy and verifiable. The fully disinfectable surface enables hygienic operation at minimal cost.

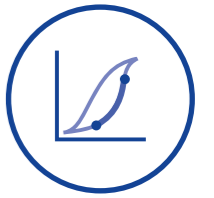

PEEPfinder®

gold standard bedside lung diagnostics

The PEEPfinder® makes the determination of the optimal PEEP range as easy as setting the respiratory rate.

It is considered an established fact that the cyclic collapse and reopening of lung areas in patients with ALI significantly damages the pulmonary tissue and that alveolar cycling of lung areas in particular represents an independent risk factor for higher mortality.

The PEEPfinder® can be used to optimise the settings of the ventilator, thus supporting lung-protective ventilation. The manoeuvre is performed in a secure window and can be combined with a preoxygenation function. The expanded, quasi-static PV tool supports the user's assessment of stress and strain. Intelligent algorithms and extensive safety features make it easier to determine the elastic properties of the lungs. A number of evaluation options are available for this purpose. Graphic evaluation support for detecting inflection points, stress indices, and storage of up to 10 reference loops facilitate the straightforward implementation of lung-protective ventilation.





Tools to assist the weaning process

there are no simple answers when weaning fails

In the majority of ventilated patients,ventilator weaning is quick and can be successfully achieved by simple strategies. However, there is a steady rise in the number of ventilated patients that cannot be weaned off the ventilator or where the weaning process is very prolonged.


40 % of all ventilated patients undergo difficult or prolonged ventilator weaning, which takes up almost 50 % of intensive

care time. Often, these are patients with severe respiratory dysfunction, where comorbidity makes the weaning process more difficult.The necessary weaning strategy is complex, demanding and allows no simple answers. In addition to special modes for simple weaning, there are numerous tools and indices available for continuously assessing the weaning process and for the standardised assessment of weaning and extubation readiness.



Weaning modes

The right choice of ventilation type has high significance in the weaning concept and influences the duration and success of weaning. In addition to the whole range of conventional ventilation modes, elisa 300 also has two special ventilation types for efficient weaning of standard ventilation patients. Spontaneous breathing activity, necessary ventilation pressure for mandatory and spontaneous breathing activities, trapping risk, and lung parameters are continuously recorded, assessed and used to adjust the ventilation parameters. Adaptive Lung Protection Ventilation (ALPV®) takes lung protective protection rules into account and guarantees the necessary CO₂ elimination.ALPV® can be maintained throughout the entire period of ventilation without changing the ventilation mode or adjusting the ventilation parameters.

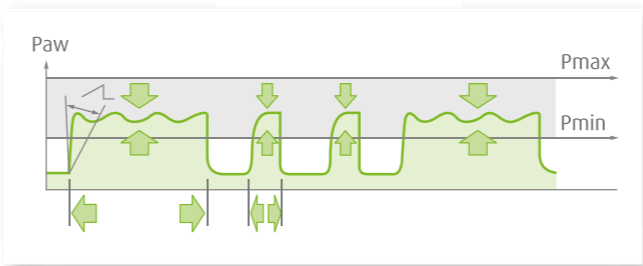


Adaptive Lung Protection Ventilation continuously adapts to the weaning situation.

Weaninganalyzer®

A huge challenge in weaning is to establish the right time for weaning readiness and extubation.The fact that up to 16 % of extubations are unplanned as so-called self-extubations with subsequent ventilation no longer being required in about 50 % of these patients illustrates the importance of the right time for planned extubation. The Weaninganalyzer® contains protocols for daily standardised determination of weaning readiness ("ready to wean") and extubation readiness ("ready to extubate").By monitoring clinical situations and assessing measurement values, daily SAT tests can be performed more easily, thus helping to reduce complications, reintubation rates, days in intensive care and treatment costs.

Fastwean®

Fastwean® allows measurement values relevant to weaning to be assessed at a glance.Whether RSBI or Negative Inspiratory Force – the measurement values are continuously displayed and assessed using a 'traffic lights' display.

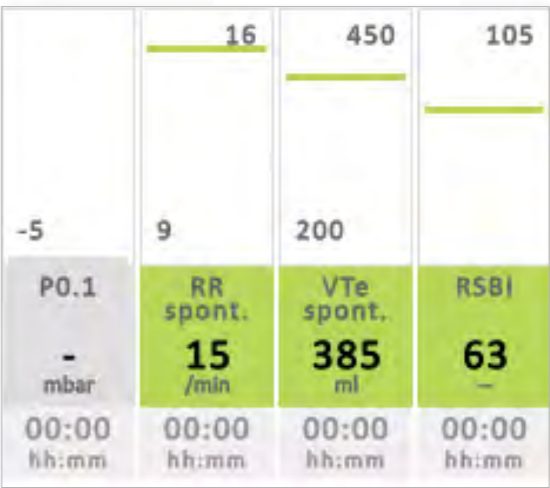

Options & choices

our modular system at a glance



Highflow O₂

High-flow oxygen therapy (HFOT) is considered a supplement to non-invasive ventilation or is used in cases where conventional oxygen therapy does not provide adequate oxygenation. It involves offering a continuous flow with individually adjusted oxygen supply via a special nasal cannula.



CPR mode

Special emergency mode for ventilation in resuscitation situations.



ALPV®

The ALPV mode combines the previous advantages of hybrid closed-loop ventilation with the current requirements of lung-protective ventilation. The pressure-controlled ventilation with volume guarantee (comparable to dynamic BiLevel) is combined with pressure-supported spontaneous breathing with volume guarantee (dynamic PSV) in such a way that a tidal volume of 6 ml/kg of ideal body weight results as the target value for mandatory and pressure-supported spontaneous breathing. At the same time, the device continuously monitors potential air trapping and offsets it as necessary. ALPV® is used as a weaning mode and generalist mode.



PAPS® Proportional Adaptive Pressure Support

In contrast to the fixed pressure support with PSV, a spontaneously breathing patient receives proportional pressure support with PAPS. The effective pressure support is based selectively on the respective increased elastic and restrictive resistance values. A special algorithm determines the current work of breathing based on elevated flow and stretch resistance in every breath and regulates the selective pressure support for compensation.



Loop package

Up to six selectable loops form the basis of differentiated assessment and derivation of treatment decisions. At the same time, up to 10 reference loops can be saved and displayed to compare with the current ventilation situation for diagnosis.



PEEPfinder®

Thanks to state-of-the-art sensor technology and its high-resolution sampling rate, the PEEPfinder® features algorithms for the reliable determination of inflection points to establish the necessary PEEP and ventilation range. The intuitive display allows a verifiable review of measuring values, transparent PEEP settings, and the assessment of stress indices as well as static compliance.



Weaninganalyzer®

The Weaninganalyzer® accurately displays the patient's weaning process and offers a reliable forecast for initiating the weaning process and extubation readiness based on daily trials and real-time data.



Mains-independent power supply

Additional batteries and an external charger allow off-grid operation for a period of at least four hours.



Automatic patient detection APD

As an additional safety function, users can activate the automatic patient detection (APD) feature on the configuration level to make it available. This prevents inadvertent switching to the standby function or turning the ventilator off as long as a patient is connected.



Hygiene function

To reduce the risk of nosocomial (hospital-acquired) infection, the ventilator's hygiene management function monitors the timely replacement of accessories that are in direct contact with the patient (tubing system, valve bar, suction system, HME filter, and nebulizer head). Monitoring and display follow the respective department requirements without the need for complex RFID chips or expensive breathing circuits.



WOBOV® Work Of Breathing Optimized Ventilation

WOBOV is a generalist mode that takes promoting spontaneous breathing, sufficient minute ventilation, an energetically optimal breathing pattern and compliance with specific lung protection rules into account. It continuously calculates the energetically optimal breathing pattern and adjusts the ventilation control (modified Otis formula) accordingly. If the ventilation is still insufficient, WOBOV gradually steps up mechanical support or the algorithm compensates the deficit up to the specified minute volume as needed.

